


Capillary Tubes and Laboratory Supplies for X-Ray Diffraction Use

Kodak Direct Exposure Film especially developed for X-ray diffraction analysis. Available in 35mm x 20m rolls and in 50 sheet boxes.

Polaroid XR-7 High Speed X-Ray Diffraction Cassettes for Laue and Precession Photography, and Polaroid 3000 Speed Film, Type 57, in single-exposure packets.

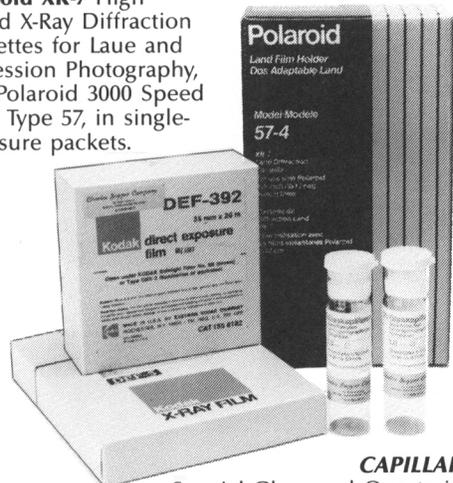

CAPILLARY TUBES
Special Glass and Quartz in 15 sizes from 0.1mm to 3.5mm O.D. Uniform 89mm length and funnel shaped at the top. 0.01mm wall thickness.

Manufacturers of Fine
X-Ray Diffraction Equip-
ment For Over 45 Years

Charles Supper Company
INCORPORATED

15 TECH CIRCLE, NATICK, MA 01760, U.S.A.
TELEPHONE: (617) 655-4610 **PD-22**

EISI Elemental and Interplanar Spacing Index

Utilizing
experimental
d-Spacings
and elemental
data...

Search
using combinations
of Elemental Data
Compound Formulae
and then ...

Identify
matches
EISILY



Now
available

for
ED
for
XRD

JCPDS - International Centre for Diffraction Data
1601 Park Lane, Swarthmore, PA 19081 U.S.A.
(215) 328-9400 Telex: 847170 Fax: (215) 328 2503 **PD-23**

Powder Diffraction File 33-34

Now available in *book form* sets.

Inorganic volume:

2934 numeric diffraction patterns of inorganic phases, metals, alloys and minerals.


Organic volume:

987 numeric diffraction patterns of organic and organometallic phases.

Crystallographic data evaluation made using NBS* AIDS83. Figure of merit for completeness and accuracy of interplanar spacings is assigned to each indexed pattern.

For further information contact:
JCPDS

International Centre for Diffraction Data
1601 Park Lane, Swarthmore, PA 19081
USA (215) 328-9400 Telex: 847170

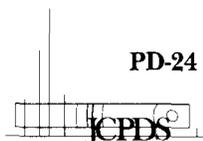

PD-24

Now available for 1989-1990...

• **Set 39 Powder Diffraction File.**

In all current media:
CD-ROM... Microfiche... Book & Mag Tape
Over 2000 new and/or updated entries in
PDF-1 and PDF-2 for:
Minerals... Metals & Alloys... Common
Phases... Inorganics & Organics

• **NIST/Sandia/ICDD Electron Dif-
fraction Data Base**

Crystallographic and chemical data on over
71,000 crystalline materials.

• **EISI Elemental and Interplanar
Spacing Index**


Rapid phase identification for electron
(& X-ray) diffractionists.

• **NIST Crystal Data File**

Crystallographic and chemical data on over
115,000 crystalline materials.

For further information contact:
JCPDS

International Centre for Diffraction Data
1601 Park Lane, Swarthmore, PA 19081
USA (215) 328-9400 Telex: 847170

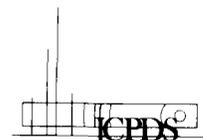

Powder Diffraction Reader Reply Card — Volume 4, Number 4, Dec. 1989

Advertiser's Information: Circle Number

PD1	PD2	PD3	PD4	PD5	PD6	PD7	PD8	PD9	PD10
PD11	PD12	PD13	PD14	PD15	PD16	PD17	PD18	PD19	PD20
PD21	PD22	PD23	PD24	PD25	PD26	PD27	PD28		

JCPDS Information: Check to Receive Information

APD Type Data Base Book Form X-ray Data Educational Packages
 Electron Diffraction Search/Match Microcomputer Search/Match
 Microfiche Package NBS Crystal Data Short Courses

Name/Title _____

Institution Company _____

Phone _____

Address _____

City _____ State _____ Zip _____

Country _____ Mail Code _____

Please enter my subscription for

Credit Cards now accepted

Powder Diffraction

Annual subscription rate: 1989

- U.S. & Canada \$39.50 (U.S. dollars)
 Library Rate \$75.00 (U.S. dollars)
 Far East Rate Write Sanyo Information System Corp., Koho Bldg. 11-8, 1-Chome,
Kayabacho Nihonbashi Chuo-ku, Tokyo 103, Japan for subscriptions for Japan, China,
Taiwan, Malaysia, the Philippines, Indonesia and Korea.

Other annual subscriptions: 1989

- \$55.00 (U.S. dollars) Library Rate \$75.00 (U.S. dollars)

Type of subscription:

- Academic Government Manufacturing Personal Other

Name/Title _____

Institution/Company _____

Address _____

Phone _____

City _____ State _____ Zip _____

Country _____ Mail Code _____

Charge: Visa MasterCard Signature _____

Number _____ Exp. Date _____

Please enter my subscription for

Credit Cards now accepted

Powder Diffraction

Annual subscription rate: 1989

- U.S. & Canada \$39.50 (U.S. dollars)
 Library Rate \$75.00 (U.S. dollars)
 Far East Rate Write Sanyo Information System Corp., Koho Bldg. 11-8, 1-Chome,
Kayabacho Nihonbashi Chuo-ku, Tokyo 103, Japan for subscriptions for Japan, China,
Taiwan, Malaysia, the Philippines, Indonesia and Korea.

Other annual subscriptions: 1989

- \$55.00 (U.S. dollars) Library Rate \$75.00 (U.S. dollars)

Type of subscription:

- Academic Government Manufacturing Personal Other

Name/Title _____

Institution/Company _____

Address _____

Phone _____

City _____ State _____ Zip _____

Country _____ Mail Code _____

Charge: Visa MasterCard Signature _____

Number _____ Exp. Date _____

TRUE MICRO DIFFRACTION

IN 3-D

Down to 11 μm Dia. or 5 μg Quantity

X-ray diffraction analysis of extremely small areas or extremely small samples can now be done with ease with the Rigaku PSPC/MDG. Three axis (ω , χ and ϕ) sample oscillation solves the problem of spotty patterns from only a few crystal grains of micro samples. The curved position sensitive detector allows simultaneous measurement of a 2θ range of 0 to 150 degrees. This minimizes the time required for data collection and produces high detection sensitivities for micro samples. A few highlights:

- 3 axis (ω , χ and ϕ) oscillation
- Simultaneous 0 to 150 degrees data collection
- Easy sample alignment
- Selection of x-ray sources
- Selection of computer processing
- Motorized X-Y-Z sliding stage (optional)
- 35 mm camera kit (optional)
- TV camera for remote viewing (optional)
- Can be added to an existing x-ray diffraction unit

Rigaku's expertise in the design and manufacture of x-ray sources and accessories is unsurpassed. For more information — give us a call !


X-ray... our only business !

Rigaku

Rigaku/USA, Inc., Northwoods Business Park, 200 Rosewood Dr., Danvers, MA 01923
Phone (508) 777-2446 ■ FAX (508) 777-3594



From now on all other powder diffractometers will be measured against this one.

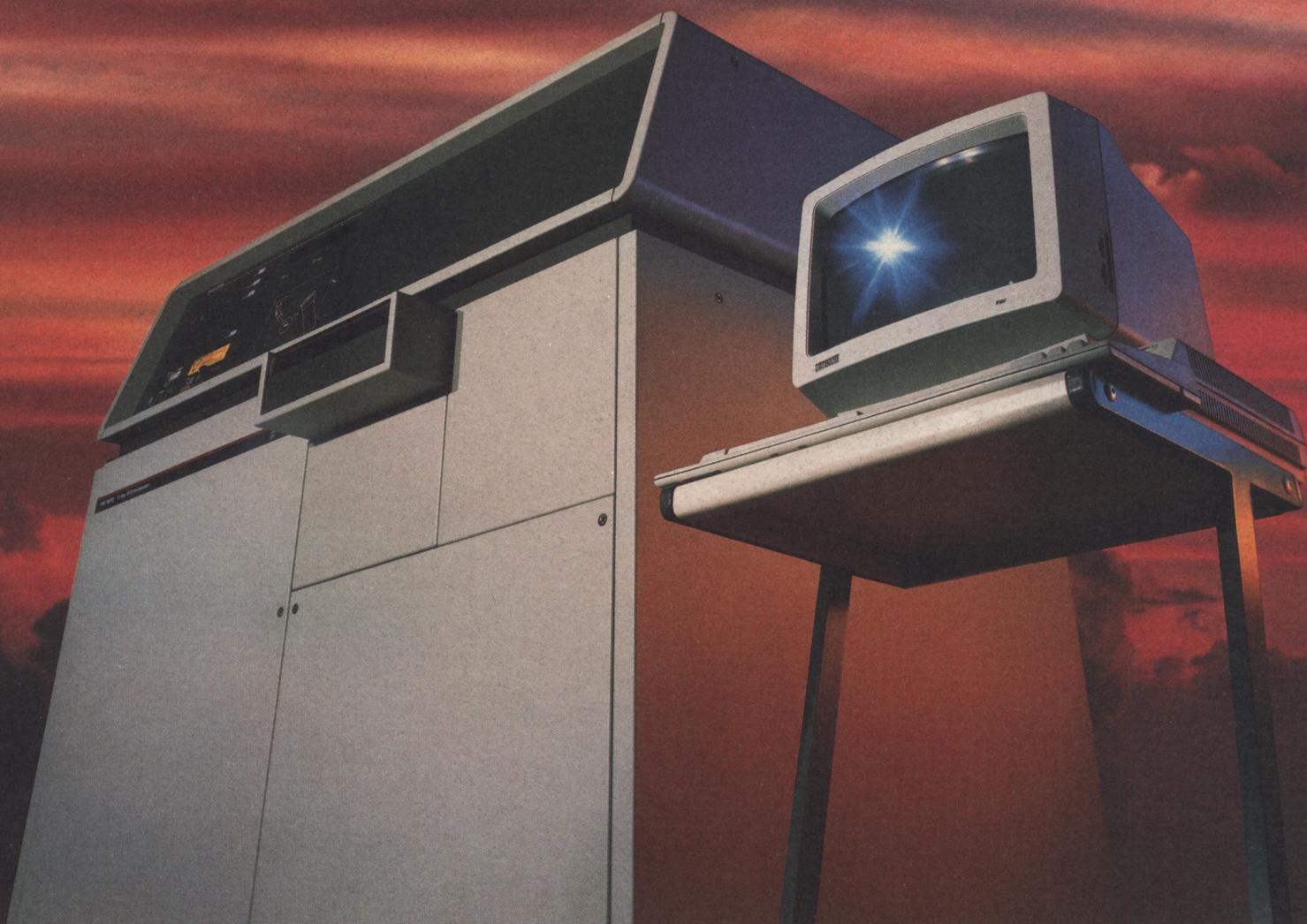

It's the new Philips PW1800.

It's the most advanced system of its type available anywhere. It's also the most versatile and the easiest to operate.

It's designed for easy, automatic measurement of a wide range of sample types in a variety of forms and quantities including powders, clays, pastes and filters. Built-in electronics control and supervise all diffractometry functions to assure reliability. What's more, the PW1800 actually monitors itself to see that it is functioning properly. And alerts the operator if it senses trouble ahead. In fact, everything about this new system has been engineered to simplify operation.

Even its appearance underscores its uniqueness. It's totally enclosed, totally integrated to assure complete radiation and contamination protection. Which means it can be installed safely and quickly in the laboratory or on the production line.

And among the things you don't see in this unique system is the newly developed, permanently-aligned goniometer that combines ease of operation with high precision and outstanding analytical performance. Built-in robotics for sample handling. And a new, highly-efficient generator that not only requires less space, but actually reduces power consumption as well.

It's the new Philips PW1800. It's the new standard of the industry. And it's backed by the worldwide support resources of the Philips organization.

For complete details, contact:

Philips Electronic Instruments, Analytical X-Ray Group,
85 McKee Drive, Mahwah, NJ 07430, Tel. (201) 529-3800.

(Outside U.S.A.) Philips Scientific & Analytical Equipment,
Industrial & Electro-acoustic Systems, Lelyweg 1,
7602 EA Almelo, The Netherlands. Tel (31) 5490-39911.



PHILIPS

PD-28